



– Cognitive Rehabilitation Software –

System requirements

- Microsoft Windows XP or later operating system
- Hardware meeting the requirements of the OS
- 1,2 GB of free hard disk space
- Display adapter and monitor (800x600 resolution and 16 bit colours minimum)
- Sound device
- Mouse or other pointing device.

Installing FORAMENRehab

Before you install:

- Check the system requirements.
- Read this instruction.

Note:

- Administrative rights are required for the installation of the FORAMENRehab software. Consult your personal support adviser if necessary.



– Cognitive Rehabilitation Software –

Detailed installation instruction

Downloading

- Download the install package from our web site: www.foramenrehab.info/index.php?page=download.htm

Starting installation


- Close all applications
- Browse to the folder where you saved the installation file
- Double click the file icon.

Language	File
English (UK)	FORAMENRehab0x0809.msi
Swedish (Finland)	FORAMENRehab0x081D.msi



– Cognitive Rehabilitation Software –

Welcome

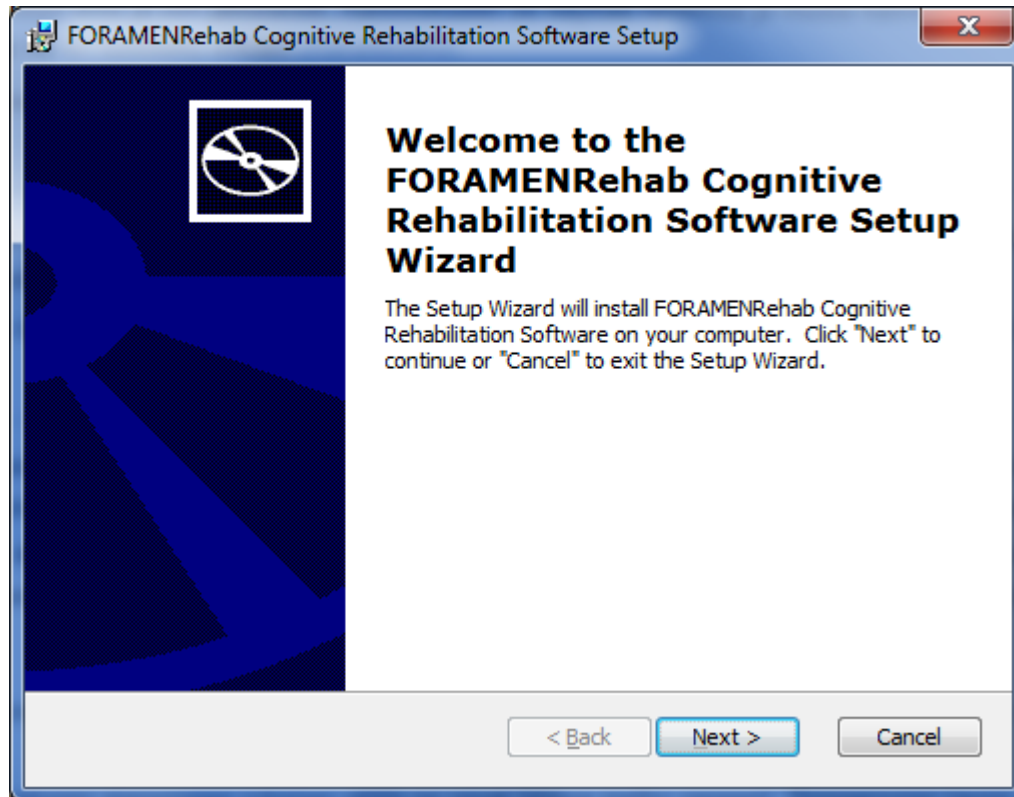


- Click the Next button.



– Cognitive Rehabilitation Software –

Select Installation Folder

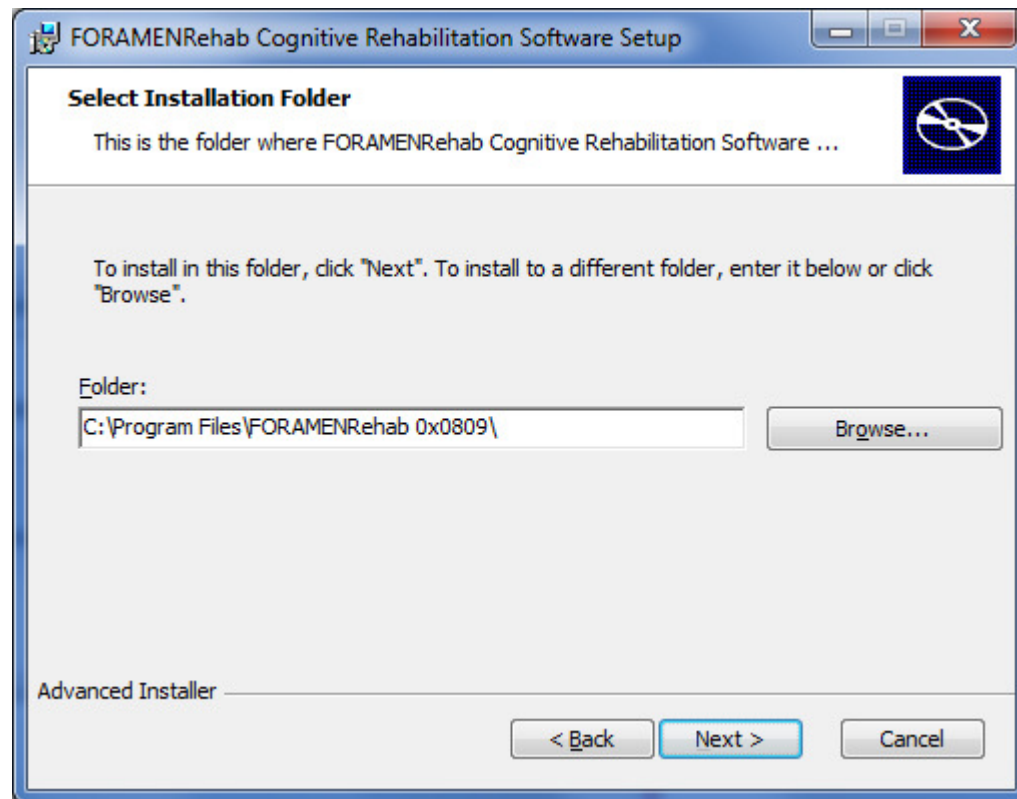


- Select the folder where FORAMENRehab will be installed (the default folder recommended)
- Click the Next button.



– Cognitive Rehabilitation Software –

Ready to Install

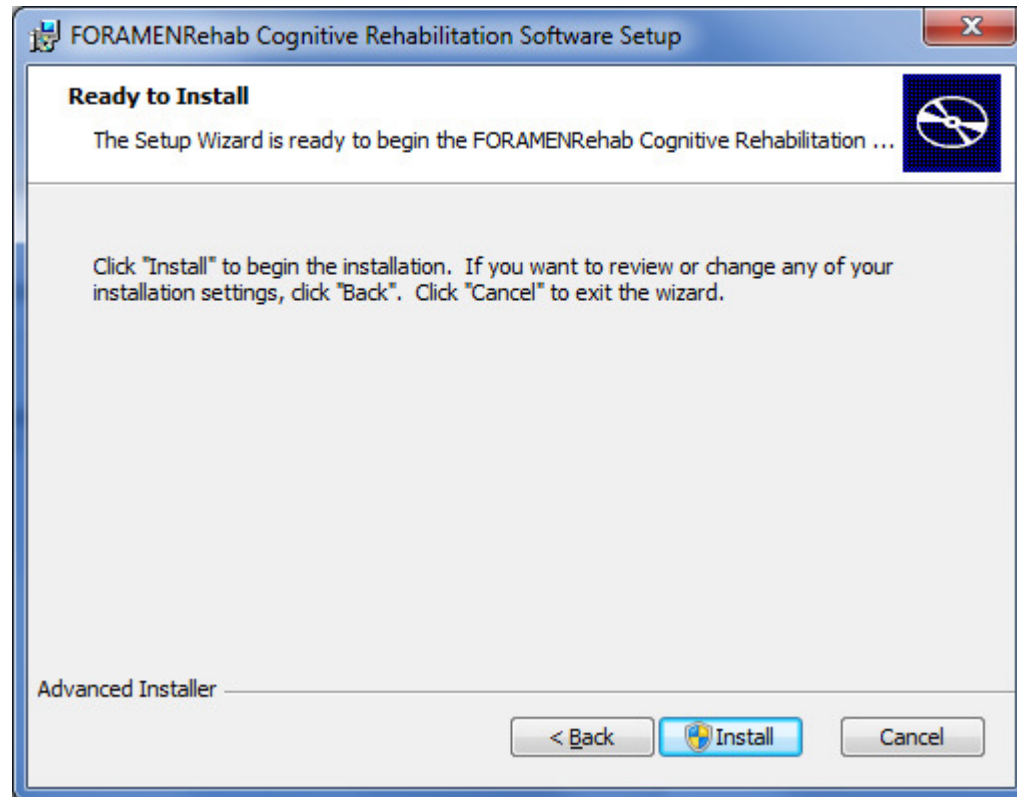


Click the Install button.



– Cognitive Rehabilitation Software –

Installing FORAMENRehab

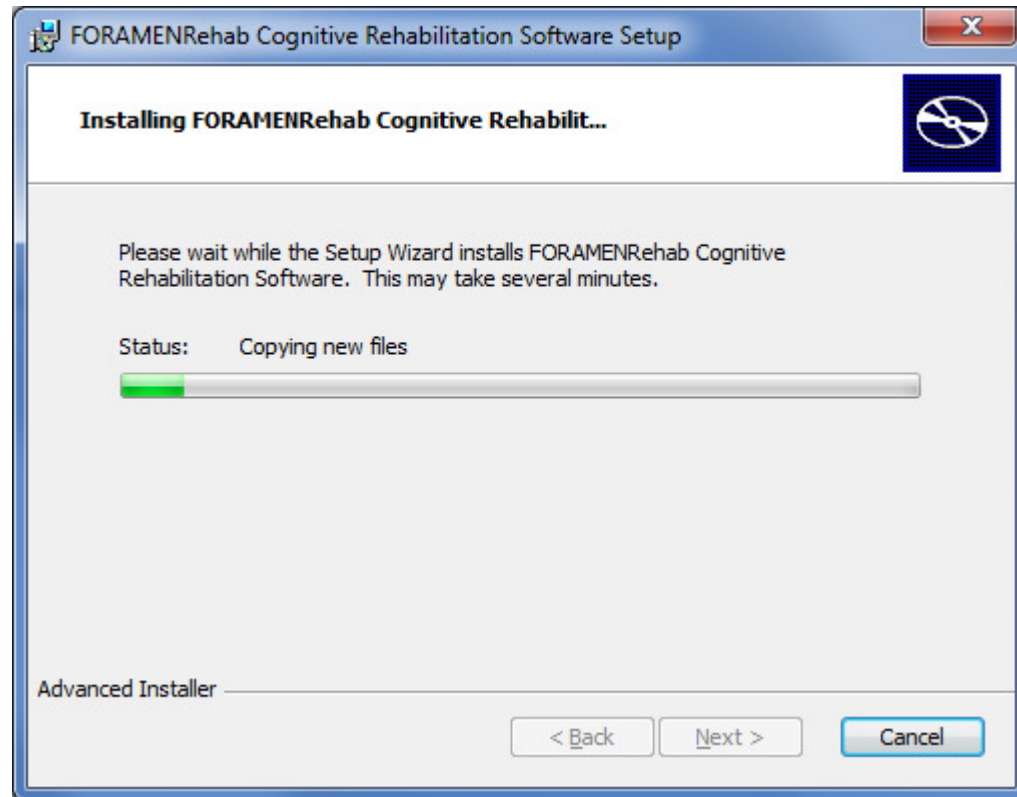


- Files are being copied into your hard disk.



– Cognitive Rehabilitation Software –

Installation complete

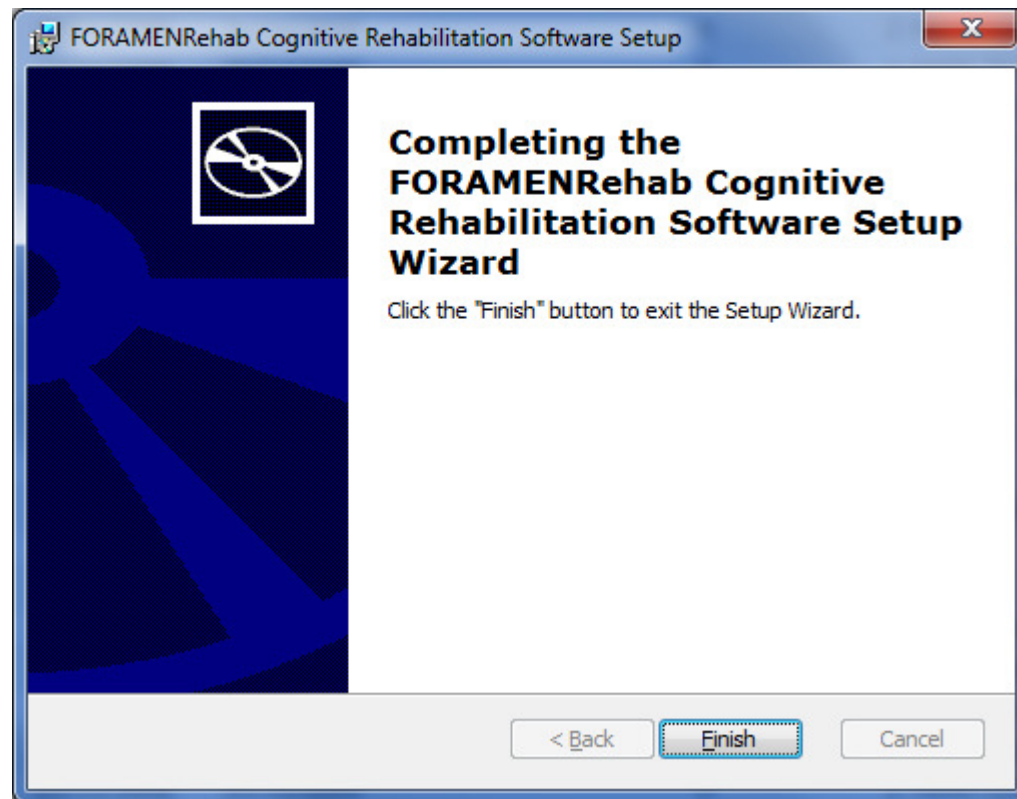

- Click the Finish button.

Grant modifying permission on the "User Files and Folders" folder to all FORAMENRehab software users of the computer:
Start-> All Programs -> FORAMENRehab Cognitive Rehabilitation Software -> Accessories -> User Files and Folders. That enables the FORAMENRehab software to save eg. results in this folder.

FORAMEN



– Cognitive Rehabilitation Software –

Activating FORAMENRehab

To run the rehabilitation task applications, the FORAMENRehab module(s) must first be activated. Administrative rights are required. Consult your personal support adviser if necessary.

One licence entitles you to use the FORAMENRehab module in one (1) workstation. Your personal/organizational licence list can be found in our web site. Click the My Licences shortcut of the FORAMENRehab Accessories group and login the Licence and Activation Service with your account.

Note:

If no Internet connection is available, you can activate off-line. Please send us the serial number found in the About window (Help -> About)

by e-mail:

support @ foramenrehab.info


The activation key and instructions are then sent to you.

To activate a FORAMENRehab module on-line, click the module-specific shortcut of the FORAMENRehab Activate group.



– Cognitive Rehabilitation Software –

Activation program starts

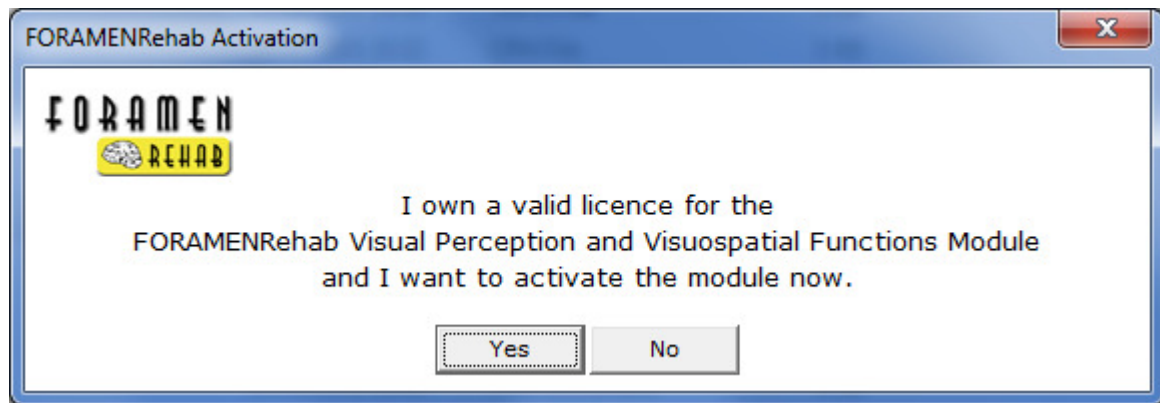


- Click the Yes button.



– Cognitive Rehabilitation Software –

Internet connection available

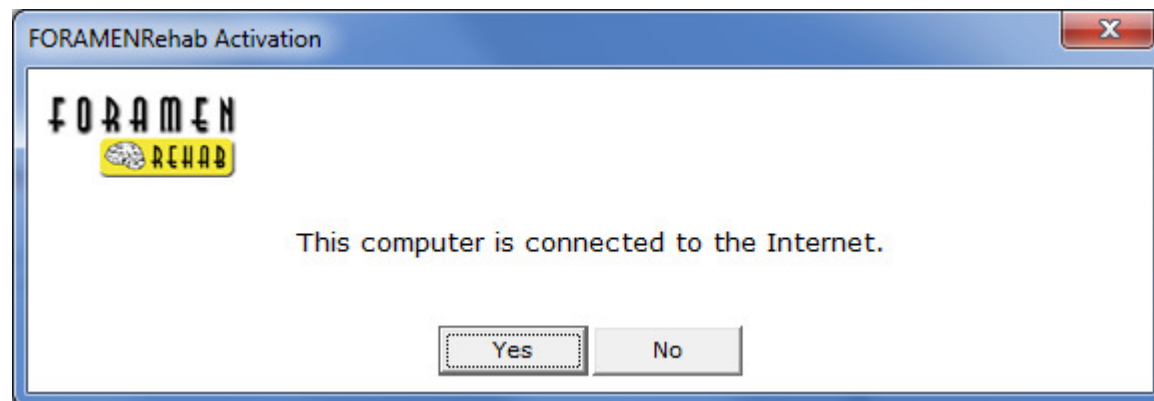

- Click the Yes button.



– Cognitive Rehabilitation Software –

The default browser is opened



- Cognitive Rehabilitation Software -

Licence and Activation Service

LOGIN

User Account:

Password:

Login

- Login the service with your username and password:

Username: [in the e-mail]


Password: [in the e-mail]

Once you have purchased a licence, a username and a password for the Licence and Activation Service are delivered to you.



- Cognitive Rehabilitation Software -


Licence list




- Cognitive Rehabilitation Software -

Licence and Activation Service


My Licenses (Visual Perception and Visuospatial Functions)

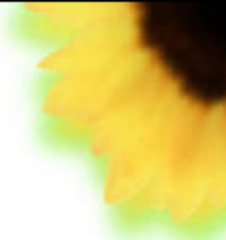

My Licenses


Visual Perception
and Visuospatial
Functions

License owner: Rehab Center, Inc.
Address: Rehab Center, Inc. 2007 Rehab Avenue Springfield, MO 65802 U.S.A.
Contact Person: Average Jones
Phone: 123-456-7890
E-mail: average.jones@rehabcenter.com


	License	Type of License	Expires	Activations Allowed	Activation Code
	9999VisualEN	1st License	N/A	1	<input type="button" value="View the existing code"/> <input type="button" value="Create a new code"/>

- Click the Create a new code button of the licence you wish to activate.



– Cognitive Rehabilitation Software –

Activation code created

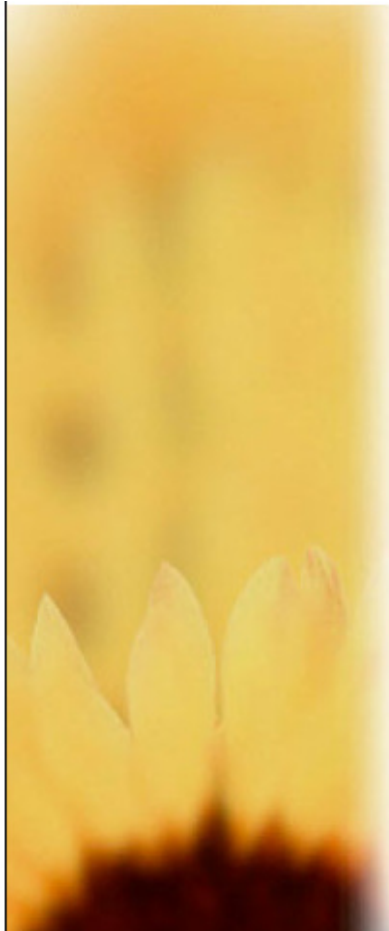

Your activation code for the selected license:

```
55612636086091850243913988912473229423966952248347717211803150707620219814004558447  
55812385213125180644106423942041744318726775473174001068693445230363229
```


Activate Now!

[Back to My Licenses](#)


- Click the Activate Now! button

A. Using Microsoft Internet Explorer:

- Click the Activate Now! button.
- The browser may ask you to allow access to the Clipboard.

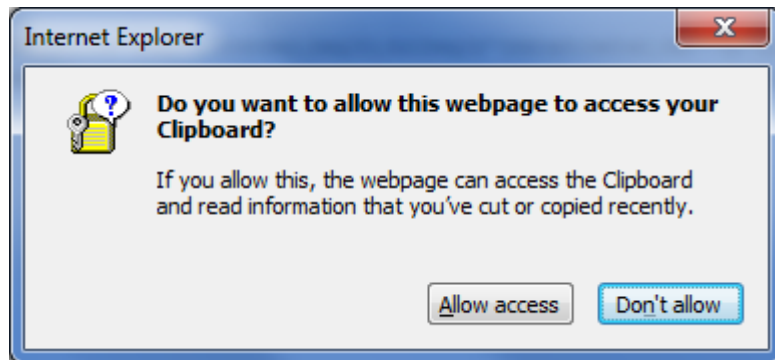

- Click the Allow access button.

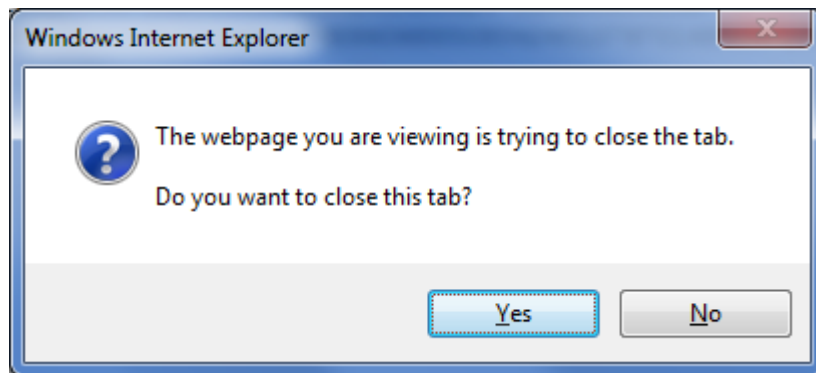




– Cognitive Rehabilitation Software –

- Click the Yes button.

B. Using other browsers:

- Type Ctrl-C.
- Close the browser.

FORAMENRehab module activated

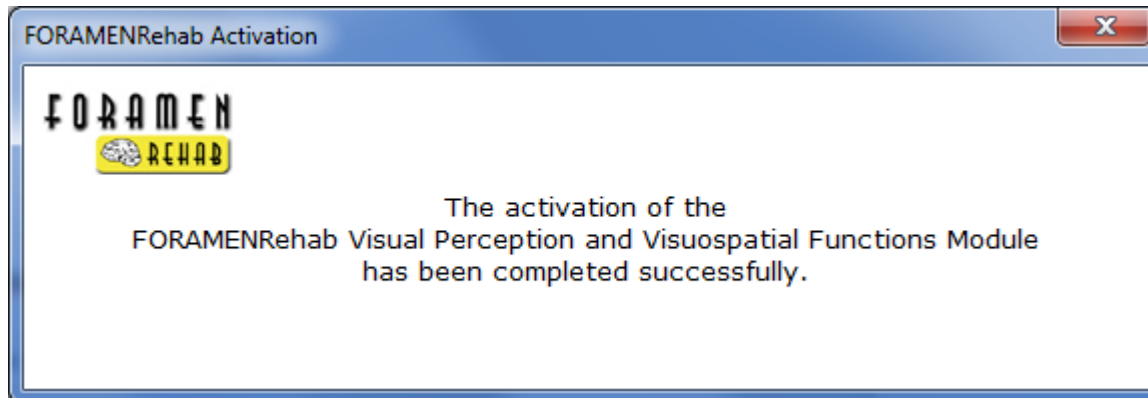

- The activation is now complete. Close the program by clicking Close Window button.



– Cognitive Rehabilitation Software –

Uninstalling FORAMENRehab


Administrative rights are required to uninstall. Consult your personal support adviser if necessary.

Uninstall the FORAMENRehab software by using the following Control Panel tools:

- Windows XP: Add/Remove Programs
- Windows Vista: Programs and Features
- Windows 7: Programs -> Uninstall a program

- Select FORAMENRehab Cognitive Rehabilitation Software
- Click the Uninstall button.

Confirming uninstall

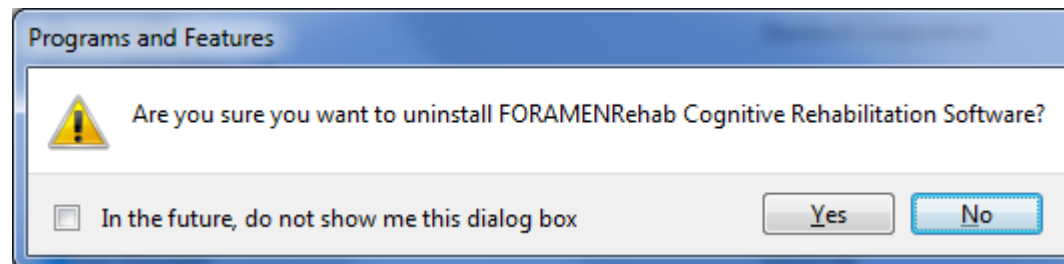


- Click the Yes button to uninstall.




– Cognitive Rehabilitation Software –

Preparing to remove

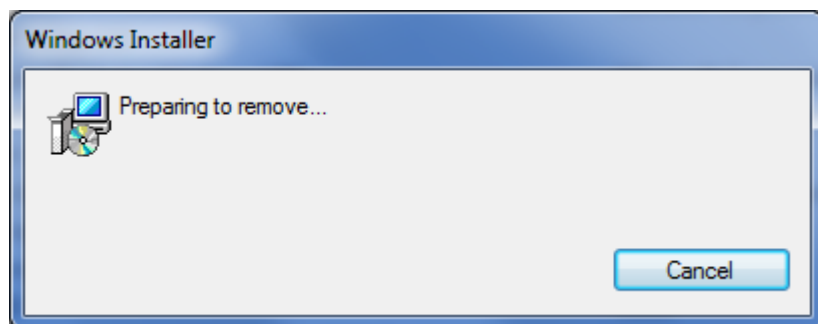

- Please wait.

Uninstalling

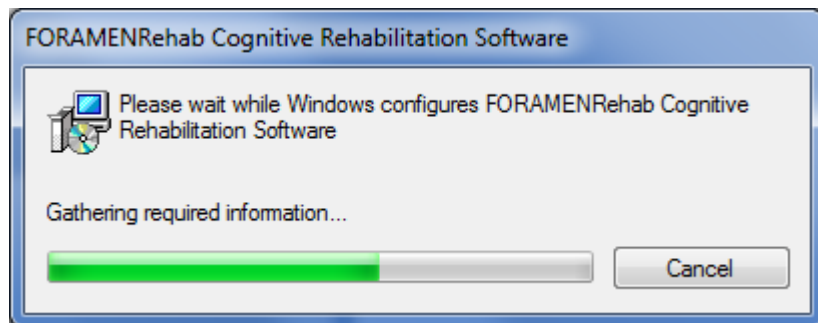



– Cognitive Rehabilitation Software –

After uninstalling

Some folders and files must be removed manually.

Depending on the OS and language version e.g.:

Language	Windows XP	Vista, Windows 7
English (UK)	C:\Documents and Settings\All Users\Documents\ FORAMENRehab 0x0809	C:\Documents and Settings\All Users\Documents\ FORAMENRehab 0x081D
Swedish (Finland)	C:\Users\Public\Documents\ FORAMENRehab 0x0809	C:\Users\Public\Documents\ FORAMENRehab 0x081D



– Cognitive Rehabilitation Software –

Quick start guide

Evaluation

Login the software with the built-in demo-account:
User account: DEMO
Password: DEMO
(NOTE! Type the password in UPPERCASE letters.)

Evaluation: see model performances, available parameters and results. The Start button is disabled until the FORAMENRehab module in question has been activated.



– Cognitive Rehabilitation Software –


Rehabilitation

To use the full-featured rehabilitation software a valid licence is required. Once you have purchased a licence, a username and a password for the Licence and Activation Service are delivered to you. Activated software can be used for rehabilitation: clicking the Start button starts the performance of a task.

Overview





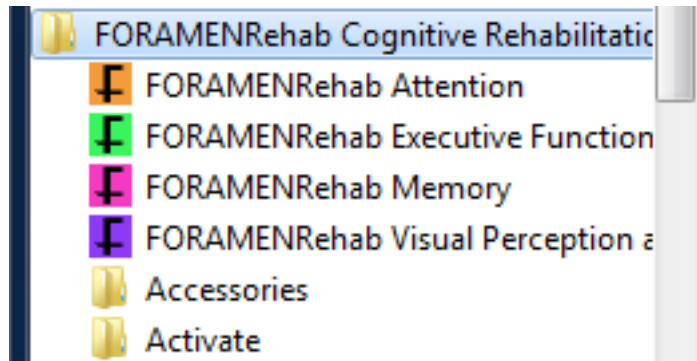

FORAMENRehab Program Groups:

I. The FORAMENRehab Program Group contains the following shortcuts:






– Cognitive Rehabilitation Software –







Shortcut	Function
 FORAMENRehab Attention	Start a module for cognitive rehabilitation of disorders of attention.
 FORAMENRehab Executive Functions and Problem Solving	Start a module for cognitive rehabilitation of disorders of executive functions and problem solving.
 FORAMENRehab Memory	Start a module for cognitive rehabilitation of disorders of memory.
 FORAMENRehab Visual Perception and Visuospatial Functions	Start a module for cognitive rehabilitation of disorders of visual perception and visuospatial functions.
 Accessories	Open a subgroup containing shortcuts to various accessories of the FORAMENRehab software.
 Activate	Open a subgroup containing shortcuts to start activation of each module separately.

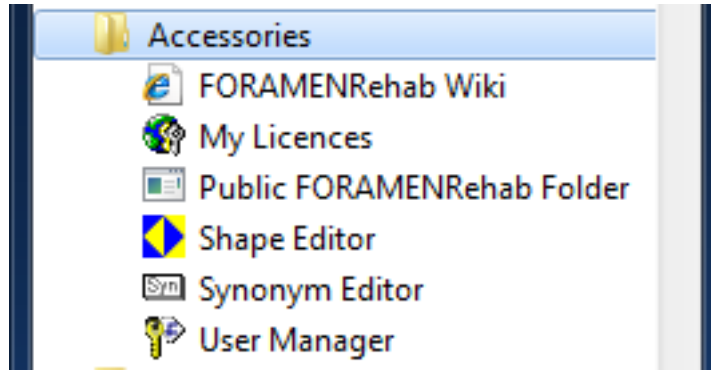

– Cognitive Rehabilitation Software –

II. The FORAMENRehab Accessories subgroup contains the following shortcuts:








Shortcut	Function
 FORAMENRehab Wiki	Open the FORAMENRehab Information Portal
 My Licences	Open the FORAMENRehab Licence and Activation Service.
 Public FORAMENRehab Folder	Open the Public FORAMENRehab Folder.
 Shape Editor	Manage the Shape Library used by construction tasks of the Visual Perception and Visuospatial Function Module.
 Synonym Editor	Manage the Synonym Dictionary used by the Memory Module.
 User Manager	Start the FORAMENRehab User manager. Password: rehab

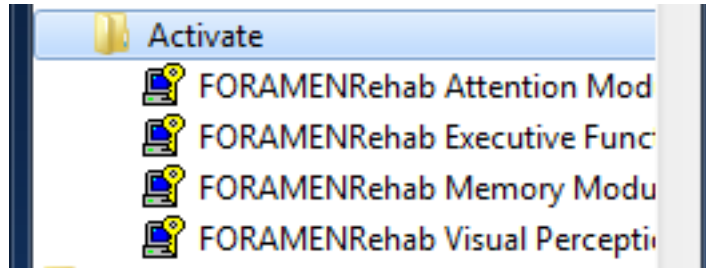

– Cognitive Rehabilitation Software –

III. The FORAMENRehab Activate subgroup contains the following shortcuts:



Shortcut	Function
 FORAMENRehab Attention	Activate the FORAMENRehab Attention Module.
 FORAMENRehab Executive Functions and Problem Solving	Activate the FORAMENRehab Executive Functions and Problem Solving Module.
 FORAMENRehab Memory	Activate the FORAMENRehab Memory Module.
 FORAMENRehab Visual Perception and Visuospatial Functions	Activate the FORAMENRehab Visual Perception and Visuospatial Functions Module.

- A valid licence is required for the activation.
- A username and a password for the FORAMENRehab Licence and Activation Service are delivered once a licence has been purchased.



– Cognitive Rehabilitation Software –

Backups

Some files and folders maintained by the FORAMENRehab software should regularly be backed up (eg. the results). These files and folders can be found in the "Public FORAMENRehab Folder".

"The Public FORAMENRehab Folder" includes:

- File: synonyms.xml (The Synonym Dictionary)
- Files: Kuviot2.dat, Kuviot3.dat, Kuviot4.dat, Kuviot5.dat, Kuviot6.dat (The Shape Library)
- Files: Module1.cph, Module2.cph, Module3.cph, Module4.cph (Activation keys for the FORAMENRehab modules; optional)
- Folder: Registry (The FORAMENRehab User Registry)
- Folder: Data (The user-specific FORAMENRehab result folders).



– Cognitive Rehabilitation Software –

The location of the "Public FORAMENRehab Folder" folder is e.g.

Language	Windows XP	Vista, Windows 7
English (UK)	C:\Documents and Settings\All Users\Documents\ FORAMENRehab 0x0809	C:\Documents and Settings\All Users\Documents\ FORAMENRehab 0x081D
Swedish (Finland)	C:\Users\Public\Documents\ FORAMENRehab 0x0809	C:\Users\Public\Documents\ FORAMENRehab 0x081D

Use the Public FORAMENRehab Folder shortcut in the FORAMENRehab Accessories group to locate the folder.



– Cognitive Rehabilitation Software –

Miscellaneous

Task applications and multi-tasking

Some FORAMENRehab task applications are time-critical. To optimize the performance, close all other Windows applications.

Regional Settings

English (UK) is recommended as the regional setting of Microsoft Windows. You may use other settings as well, but settings of an individual user should not be changed once FORAMENRehab has been used for the first time. This ensures the readability of the results.

Credits

The software uses images of Corel GALLERY 1 000 000 -library, copyright owned by Corel Corporation. You may not reproduce or distribute the images in any way.